
11.15.10

L01-429 Rev. A

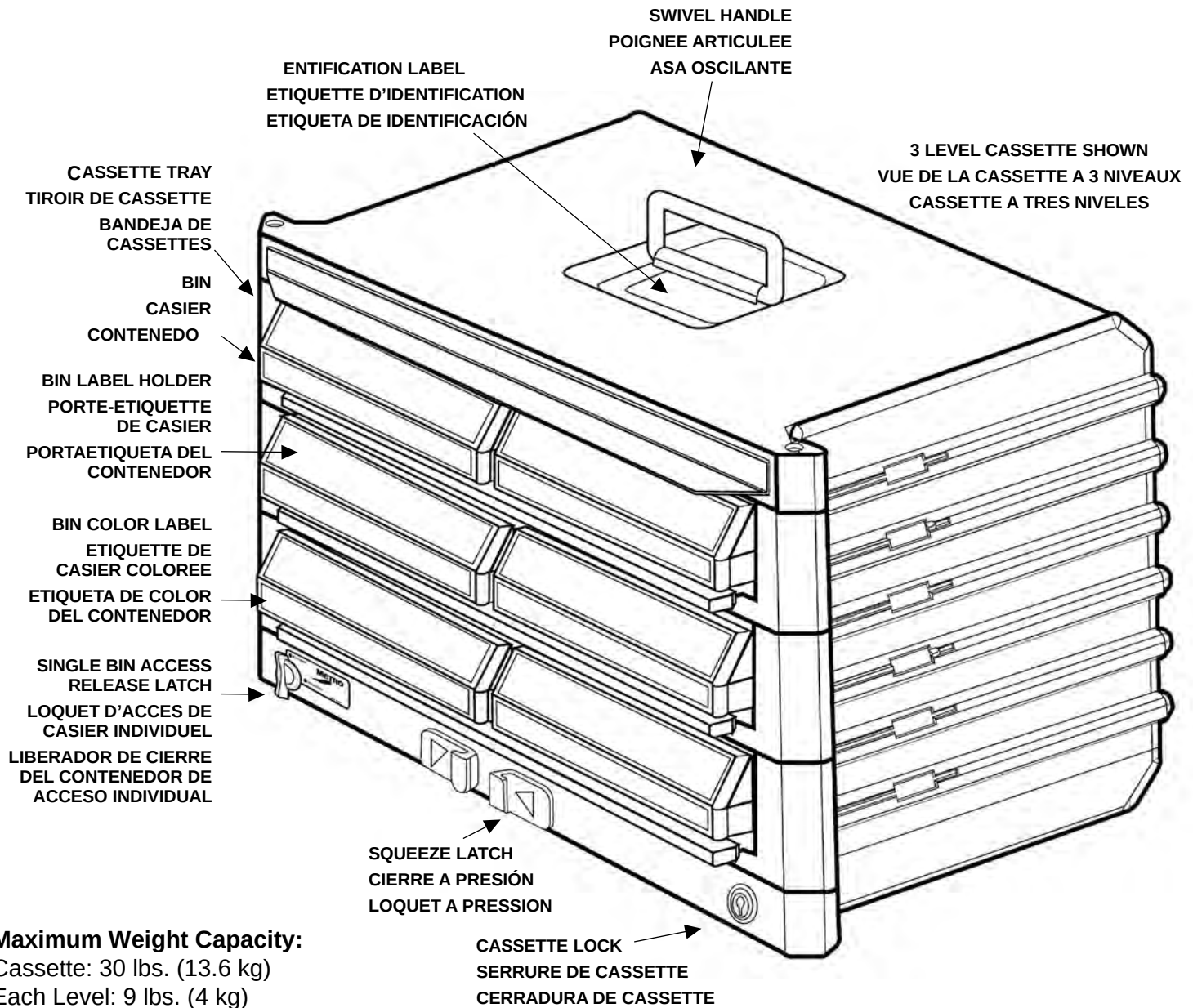Maximum Weight Capacity:

Cassette: 30 lbs. (13.6 kg)
Each Level: 9 lbs. (4 kg)

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Cassette: 13.6 kg
Niveau individuel: 4 kg

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Cassette: 30 lbs (13,6 kg)
Cada nivel: 9 lbs (4 kg)


Replacement Parts on Page 3
Cleaning Instructions on Page 4

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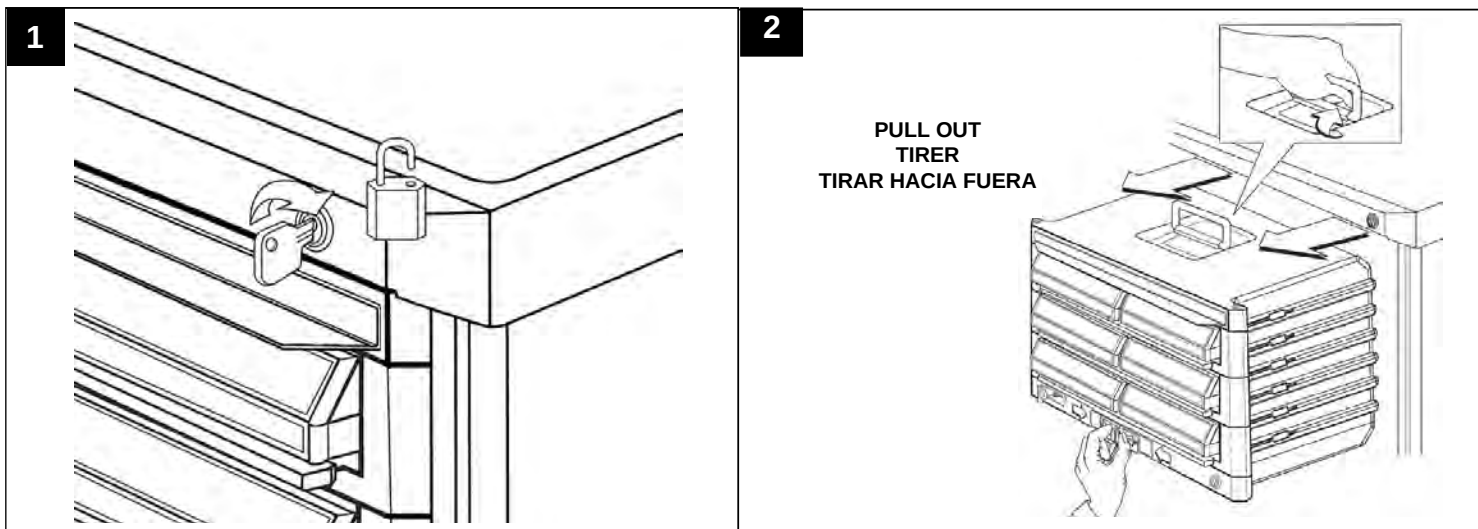
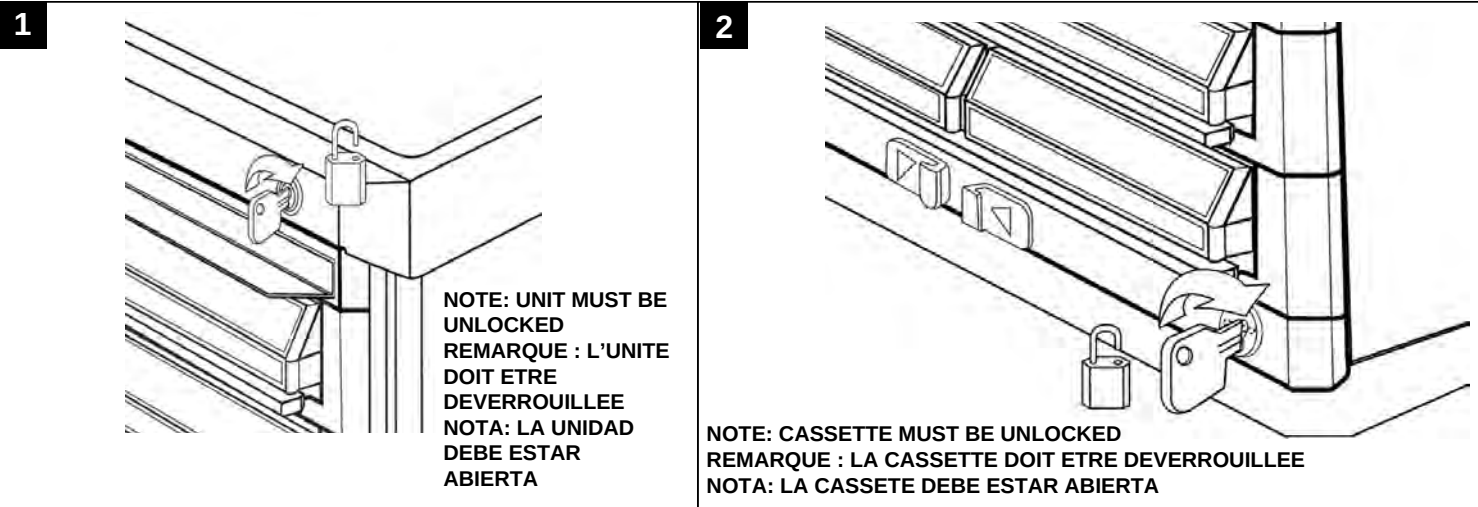
FULL LEVEL ACCESS
ACCESS RANGEE COMPLETE
ACCESO AL NIVEL COMPLETO

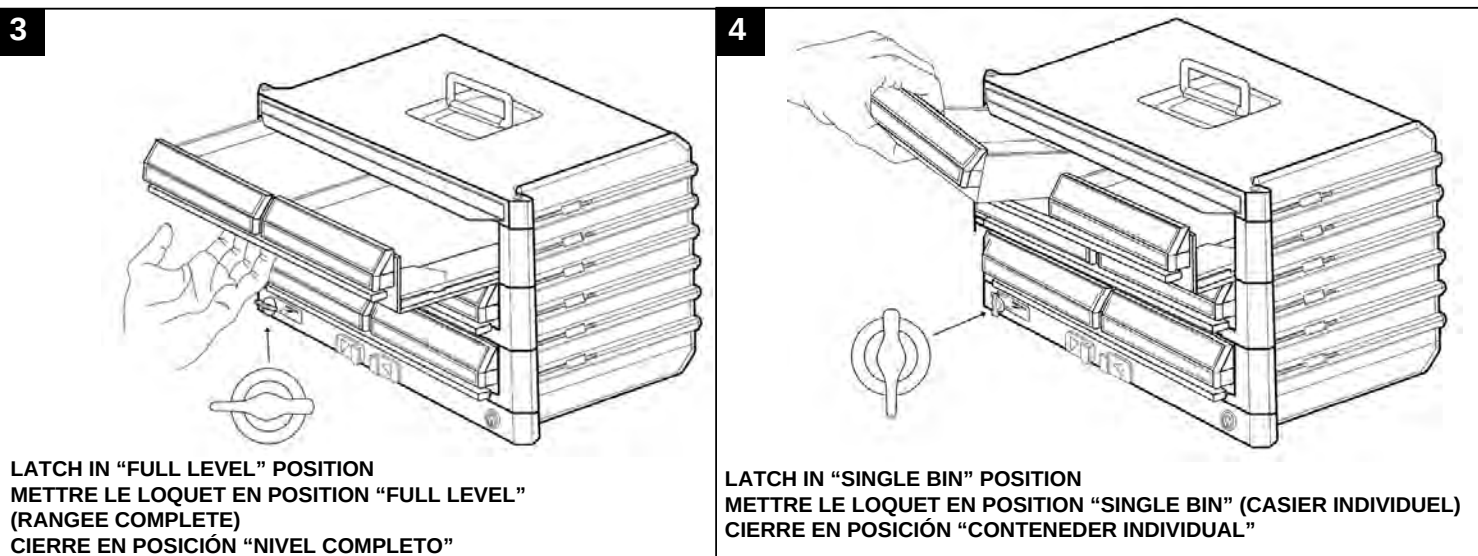Cassette Body


Item	Description	Qty.	Part No.
A	Cassette ID Label	50 Pack	RPCASID-50
B	Latch Spring	1	RPC10-016
C	Latch with Screw	1	RPCAS-LTCH
D	Lock & Cam Assy.	1	RPCAS-LOCK
E	SBA Latch & Cam	1	RPCAS-SBALK
L	Cassette Tray	1	RPSXR-1020

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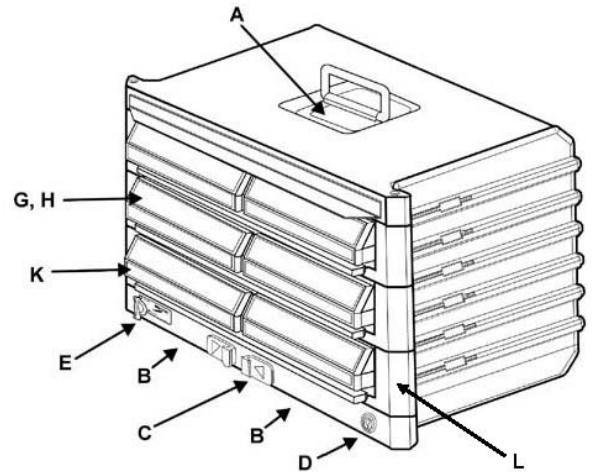
* Cassette Label Kit includes: (18) 3", 4 1/2", 6" & 8" Bin Cards; 2 Sheets bin labels (1-18/sheet 3", 4 1/2" & 6" bins; 1-33/sheet for 8" bins) in each color & 4 Cassette ID Labels.

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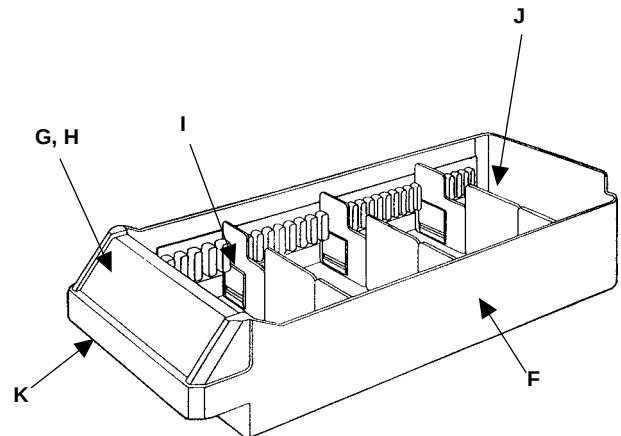CLEANING INSTRUCTIONS

Starsys cart bodies and some accessories are constructed of durable advance composite of polymers. A soft cloth, mild soap and water solution can be used to clean lightly soiled surfaces. After cleaning, wipe dry with a clean soft cloth.

For heavily soiled areas, a soft brush and solvent or emulsion-based cleaners may be employed. Always insure that the cleaner is recommended for use on plastics and follow any special instructions from the manufacturer.

The use of abrasive cleaners is **NOT** recommended.

Please read the instructions for recommendations on correct cleaning materials and procedures.

RECOMMENDED TYPES OF CLEANING PRODUCTS:

NOTE: This list is by no means complete. The listed products are only ones that have been proven effective in cleaning the cart body, drawers and accessories. It is not our intention to promote particular brands.

CLEANER BASE	DESCRIPTION	PRODUCT EXAMPLES
Phenol	Liquid detergent, disinfectant, deodorizer	Amphyl
Chlorine	Liquid disinfectant	Bleach
Alkaline	Liquid detergent Liquid detergent	Conquer II Prime
Alcohol	Antiseptic sterilizer	Isopropyl Alcohol
Iodophor	Liquid germicide detergent	Accord
Quaternary Ammonia Compound	Liquid detergent, disinfectant Liquid detergent, disinfectant Liquid detergent, germicide, deodorizer Foam cleaner, disinfectant Liquid detergent, disinfectant	FSD-4 Spartec Micro-Quat Spra-Temp Lysol 256

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